




-  Fast & accurate way to fill capsules without spillage.
-  Use Tamper for higher fill weights.
-  Stainless Steel (highly polished) for high GMP and easy cleanup.
-  Lightweight for use on analytical balances.
-  Available in Sizes 000-5 including el sizes; for rats 9, 9el; and for veterinary use 11 & 12/13.

*Free Capsules offer for sizes 000-5 excluding el & veterinary sizes. Not offered with purchase of ProFunnel Kit.

ProFunnel Kit



One funnel of each size 000-5

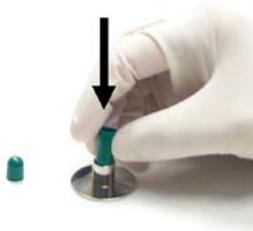

1. Separate a single capsule. Place the capsule body (Longer, narrower part of the Capsule) in the holder.




2. Place the funnel over the capsule body in the holder.



3. Dose a pre-measured amount of powder into the capsule through the funnel, **or** put the assembly on a precision balance, tare it, and dose the powder into the capsule through the funnel.



4. When the required quantity of powder has been put into the funnel, use the tamper to compress filling material into capsule body.

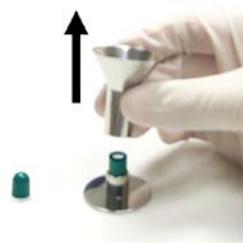

5. Slowly remove the funnel.



6. Place the capsule's cap onto the body of the capsule, pressing firmly to lock. Remove the filled capsule from the base.

Do you need to fill more than a few capsules?

ProFiller 1100



Hand-held filler up to 3200 capsules/hour. Sizes 000-5, 00el, 0el. For sizes AAA-D ask about ProFill OE.

ProFiller 3700/3800



300-hole Bench-top filler 6000 - 9000 capsules/hour. Optional vibrator & European style Orienter.

ProCoater



Enteric / Film Coat Capsule & Caplets. Enteric coat for lower intestine drug release. Film coat for taste masking or appearance.

ProMixer



Powder Mixer/Blender for uniform powder mixtures. Choose from 6 different Inox (SS 316) shells for maximum flexibility.